






### Aesthetic Product List

	Type of machine	Description
<p><b>ULTRAFORMER<sup>®</sup> III</b></p> 	<p><b>MMFU</b> Micro and Macro Focused Ultrasound</p> <p><b>MMFU</b> Micro thermal effect High peak power with Less pain by dual engine</p> <p><b>Multi-functional cartridges</b> 1.5mm cartridge (L7-1.5) 3mm cartridge (L7-3.0) 4.5mm cartridge (L4-4.5) 2mm cartridge (MF2) 6mm cartridge (MF6) 9mm cartridge (MF9) 13mm cartridge (MF13)</p>	<p><b>ULTRAFORMER III</b> is a device based on <b>MMFU</b> (Micro &amp; Macro Focused Ultrasound) technology for skin lifting and body contouring. Non-invasive UltraformerIII treatment effectively works on both face and body with no downtime. UltraformerIII has been certified it's safety in many countries and used in all over the world.</p> <p>Multi-functional cartridges assisted by MMFU (Micro &amp; Macro Focused Ultrasound) technology are specially designed for body treatment to combine with previous 1.5, 3.0, and 4.5mm cartridges in order to enhance its treatment effect. Current 4&amp;7Mhz transducers create Micro coagulation into dermis and SMAS while 2Mh transducers deliver Macro energy into subcutaneous fat to selectively remove fat cells. This innovative transducers that generates 8 times higher ultrasound energy maximize customer's satisfaction and reduce treatment time.</p> <p><b>ULTRAFORMER III</b> has various cartridges in order to treat any parts of skin and body such as face, neck, abdomen, buttocks, upper arms, knees and so on. The practitioner can select a cartridge according to the thickness of the target area. During <b>ULTRAFORMER III</b> body treatment, the membranes of the fat cells are damaged. The Fat cell content, primarily comprised of triglycerides, is dispersed into the fluid cell and then transported through the vascular and lymphatic systems to the liver. The liver makes no distinction between fat coming from <b>ULTRAFORMER</b> treatment and fat originating from consumed food. Both are processed via body's natural mechanism.</p>
<p><b>CLATUU α</b></p> 	<p><b>CLATUU ALPHA</b> Cryolipo therapy Non - Invasive Natural Body Shaping And Contouring</p> <p><b>VARIOUS APPLICATORS FOR ANY SPOT</b></p> <p>User Friends Operation System Effective Result</p>	<p><b>The Clatuu Alpha</b> is a medical aesthetic device that provides non surgical fat reduction treatment using upgraded 360° surround cooling technology.</p> <p>The Clatuu Alpha has been upgraded with more efficient colling technology, allowing an 18.1% increase in efficiency. The Clatuu Alpha is also equipped with an expandable range of smart applicators that can automatically computer recommended parameters for each treatment applicator. The upgraded system also allows users tpo adjust cooling and suction levels to provide treatments customized to the needs of each patient.</p> <p>As the Clatuu Alpha treatment induces apoptosis, it takes a little bit of time to see the result. Therefore, with Clatuu Alpha treatment you will start seeing results after 4-12 weeks.</p> <p>According to the study, normally one session (applying the 1 applicator) reduces size by 1 inch. As Clatuu Alpha has two applicators,patients reduces the size 2 inches at a time.</p>



### Aesthetic Product List

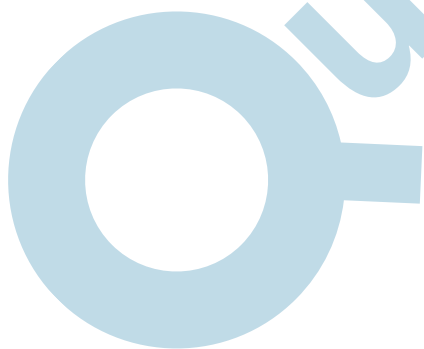
	Type of machine	Description
	<p><b>CLATUU</b></p> <p>Cryolipo therapy Non-Invasive Natural Body Shaping and Contouring Dual Applicators User Friendly Operation System Effective Result</p>	<p><b>Clatuu</b> is the non-invasive cooling of fat cells to induce lipolysis without damage to other tissue types. CLATUU fat freezing away is a safer alternative to liposuction. The purpose of CLATUU treatment is to remodel the shape of the body.</p> <p>CLATUU treatment, cryolipo therapy works on the basis that fat cells are more vulnerable to cold. That is why CLATUU treatment makes fat cell freeze selectively without damage of other tissue. The destroyed fat cell is, then, naturally released from body through metabolic process over a period of 6 to 12weeks.You might see the redness the area you get treated. It disappears within one hour.</p> <p>After the treatment, the cooled fat cells begin a process called "apoptosis" and begin to shrink. These damaged fat cells are slowly metabolized over 6 to 12 weeks by our bodies.The initial feeling when the applicator is applied is that the skin is being pulled intensely into a vacuum. The cooling process begins which might sting or feel a cold sensation for the first 10 minutes. Afterward, you feel very comfortable until the treatment is over.</p> <p>As the CLATUU treatment induces apoptosis, it takes a little bit of time to see the result. Therefore, with CLATUU treatment you will start seeing results after 4~12weeks.</p> <p>According to the study, normally one session (applying the 1 applicator) reduces size by 1 inch. As CLATUU has two applicators, patients reduce the size 2 inches at a time.</p>
	<p><b>ULFIT</b></p> <p>Micro Focused Ultrasound and Macro Focused Circular Ultrasound</p> <p><b>Multi-functional cartridges</b></p> <p>1.5mm cartridge (L7-1.5) 3mm cartridge (L7-3.0) 4.5mm cartridge (L4-4.5) 6mm cartridge (MF6) 13mm cartridge (MF13)</p>	<p><b>The ULFIT</b> is completely non-invasive, face lifting and tightening, body tightening and contouring device that provides customized treatment procedures applicable to the face, neck, and multiple areas of the body.</p> <p>Utilizing Focused Ultrasound technology with Micro Focused and Macro Focused Circular Cartridges, the ULFIT focused both uniform and perpetual ultrasound energy that induces thermal coagulation in the face and contraction in the body, excluding direct contact with or damage to the epidermis.</p>

**Aesthetic Product List**



	Type of machine	Description
 <p><b>SCIZER</b></p>	<p><b>SCIZER</b></p> <p>2 MHz - 13.0mm Scientific Solution For Size Reduction Reduce your size 2X Hand Pieces Larger Spot Size Shorter Procedure Times The Results Less pain</p>	<p>The Scizer treatment is designed to manipulate the adipocytes in the subcutaneous fat layer utilizing a non-invasive therapeutic focused ultrasound procedure. Based on the process of coagulative necrosis development, the treatment targets and reduces localized fat deposits, to ultimately and effectively improve the contours of the body.</p> <p>The Scizer, equipped with progressive uniform scanning technology, distributes an absolute and homogeneous beam profile to the Macro area of 46mm x 46mm.</p> <p>Full coverage of the emission area is achieved with a precise and equal distribution of ultrasound energy to the treatment area. The design of the projected output is symmetrical in size and power, to accomplish a higher yield of impact and results. Immediate Results and Gradual Improvement</p> <p>Immediate Results from Scizer treatment is evident immediately after procedures and continue to reflect positive development gradually over periods of time. By delivering ultrasound energy at 2MHz and 13mm deep, without the distribution of heat diffusion to surrounding tissues, the targeted fat cells are destroyed with the utmost care. The Scizer is recognized as a safe and effective method of fat reduction, and provides all users the highest level of satisfaction for treatment procedures and outcomes.</p>
 <p><b>AQUAPURE</b></p>	<p><b>AQUAPURE</b></p> <p>Aquapeel : Skin Exfoliation &amp; Moisture</p> <p>Electroporation : Trandermal Delivery</p> <p>Micro-Current : Facial Lifting</p> <p>Cooling : Skin Relaxation</p>	<p><b>Aquapeel</b> : Three different types of optimized solutions supply nourishing properties to the skin, dissolving impurities and helping consumers to achieve a healthier reflection.</p> <p><b>Electroporation</b> : Enhances transdermal delivery of positive ingredients deep within the skin by providing an effective and direct route to the inner layers of the skin.</p> <p><b>Micro-Current</b> : Muscle stimulation by micro-current therapy helps to achieve Elastin and collagen regeneration to establish a lifting effect by directing low electrical currents to the skin and facial muscles.</p> <p><b>Cooling / Heating</b> : Applying cooling at -5°C or moderate heating techniques post treatment is a complementary finishing step which actively relieves the skin of hyper sensitivity by influencing and stimulating blood circulation.</p>

**Aesthetic Product List**

	Type of machine	Description
	<p><b>REFIT</b></p> <p>Multi-Polar RF Skin Rejuvenation &amp; Tightening</p> <p>LED Custom selected wavelength of 640nm/460nm works</p>	<p><b>The Refit</b> utilizes radio frequency energy to effectively tighten and rejuvenate the skin by manipulating the tissue's resistance within the various layers of the skin. The RF energy is transformed into the thermal energy to trigger collagen development and to achieve skin rejuvenation.</p> <p>The Refit, the latest in radio Frequency technology represents the synergy between Multi-Polar RF, and LED Functions to competently achieve collagen contraction and production, facial lifting, tightening, cellulite reduction, wrinkle reduction, and skin rejuvenation treatments.</p> <p>RF technology since its introduction in early 2000, has been mastered to expose your potential with zero pain, no downtime, and entirely at your convenience</p>



**Aesthetic Product List**

	Type of machine	Description
<p><b>ULTRAFORMER</b> .....<b>MPT</b>.....</p> 	<p><b>NEW !!</b></p> <p><b>Ultraformer MPT</b> Faster treatment with new features With various coagulation form and depths enabling multimodal approaches</p> <p><b>Facial cartridges</b> 1.5mm cartridge ( 7 MHz ) 2mm cartridge ( 5.5 MHz ) 3mm cartridge ( 7MHz ) 4.5mm cartridge (4 MHz )</p> <p><b>Body cartridges</b> 6mm cartridge ( 2 MHz ) 9mm cartridge ( 2 MHz ) 13mm cartridge ( 2 MHz )</p> <p><b>Ultra Booster cartridges</b> 1.5mm cartridge ( 7 MHz ) 3mm cartridge ( 7 MHz ) 4.5mm cartridge ( 4 MHz )</p>	<p><b>Ultraformer MPT</b> Creating accurate and consistent TCP quality determines the performance of MMFU equipment. Exellent TCP quality A new ULTRAFORMER series, MPT , is safer and maximizes MMFU efficiency by divide the TCP generating method into micro-unit beyond the technology to create precise TCP quality</p> <p><b>Features :</b></p> <ul style="list-style-type: none"> <li>- Micro-Pulsed Technology Delivery of superimposed energy to the treatment area with 25 times finer TCP.</li> <li>- Multi select MMFU x 4 Normal, Micro Pulse, Circular, Micro Circular</li> <li>- Ultra Booster Handpiece Transducer with optimal minimization.</li> <li>- Infinite Expendables Infinite Expendables with 3 handpieces &amp; 10 interchangeable cartridges.</li> <li>- Reduced Treatment Time More than 2.5x faster transducers Elimination transducer return time.</li> </ul>

### Aesthetic Product List

	Type of machine	Description
	<p>LDM Ultrasound</p>	<p><b>LDM - MED</b> is a further development of the LDM appliance range from Wellcomet and is an ultrasound machine specially developed for the medical sector combining two special technologies (LDM Plus and the 1.0 MHz frequency).</p> <p>For progressive, holistic skin treatment, the LDM - MED machine offers treatment opportunities for various skin diseases and programs for non-invasive aesthetic applications. With up to 19 program settings divided into aesthetic and dermatological menu items, the LDM-MED covers the important indications, providing an effective equipment-based treatment system for use in dermatology, cosmetic surgery and general medicine.</p> <p>For the body, the LDM-MED has programs for different body applications such as cellulite, stretch marks and treatments before and after cosmetic surgery.</p> <p>In addition, the LDM-MED machine can be fitted with a special upgrade for body treatments (see Applications). Treatment is based on a specially developed body concept using the TWL method (Thixotropic Washout Lipolysis).</p>
<p>LTRA</p>	Type of machine	Description
	<p>MIIN Q-Switched Nd:YAG Laser Multi-Sub Pulsed Nd: YAG Total Solution for Tattoo &amp; Pigment Scar, Acne Scar, Rejuvenation</p>	<p><b>MIIN</b> applies the Nd:YAG Laser source and is composed of main body part, and arm part mainly. The main body part is consisted of laser module from which beam is irradiated, power module, control module and touch screen.</p> <p><b>Specialized to treat Melasma</b> MIIN Mode -4 - 8 multipulse Mode Maximum Energy -3J (ex, 500m J 6 pulses) It is treated as high energy as 3J, but it is split into multi-pulse with nanosecond. A specialized Melasma treatment program to reduce side effects.</p> <p><b>Features :</b></p> <ul style="list-style-type: none"> <li>- Uniform Beam Profile</li> <li>- DuExcellent Durability</li> <li>- Low Maintenance Cost</li> <li>- User Friendly Interface</li> </ul>